



Questions? Call 1-800-848-4670


Model #: \_\_\_\_\_  
Serial #: \_\_\_\_\_  
Date: \_\_\_\_\_

Provider Name: \_\_\_\_\_  
Phone #: \_\_\_\_\_  
Contact Person: \_\_\_\_\_

### Setting up your equipment:

- Home Use:
  - Fill bag and hang on IV Pole. (Picture 1)
  - Attach tubing to pump. (Picture 2)
  - Pump should be set to OFF before plugging into electrical outlet.
  - Plug in Input 120V, 60Hz, 200 mA
- Ambulatory:
  - Check for battery charge.
  - Fill bag. (Picture 3)
  - Attach to pump.
  - Place into ambulatory bag.

Picture 1

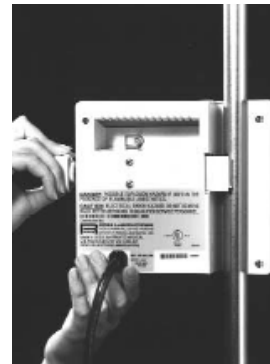

Picture 2



Picture 3



Picture 4



### Using your equipment:

- Turn pump dial to the Set Rate position, and set designated food rate (ml/hr). (Picture 4)
  - Locate correct dose and set to instructions
  - Press/Turn/Set to run.
- \* If product is not filled at bedside, it is advisable to close slide clamp before transporting.  
\* Make sure slide clamp is open for set up.

Picture 5



Picture 6



**Frequent replacement items:**

- Tubing (Picture 5)
- Feeding Bags (Picture 6)

**Special Procedures:**

- Ambulatory Pump:
  - Must have correct ambulatory bag for pump & tubing bags to work properly. Brands may differ. Check instructions for placement of pump and bag.
  - A new pump should be plugged into an A/C outlet before initial use for directed amount of time to assure full battery charge.

**Safety Issues:**

- Alarm Sounding:
  - Check for kink(s) in the tubing.
  - Check for empty bag (do not operate pump with empty bag)
  - Check for low battery
  - Check for clogging in the line, slide the clamp to assure it's completely open.
  - Check drip chamber to be sure it is not more than  $\frac{1}{4}$  full or coated with product.
  - Be sure all bags and tubing are compatible with the pump you are using.

**Maintenance:**

- The outside of the pump can be cleaned with warm soapy water.
- The pump cavity can be cleaned with a cotton swab or soft cloth and warm soapy water. DO NOT CLEAN WITH ALCOHOL.
- Pump should be turned off or recharge battery (if available) when not in use.
- DO NOT SUBMERGE OR AUTOCLAVE the pump or charger.

**Please contact your equipment provider for any questions, adjustments, or repairs.**

~

©Copyright The VGM Group 2002